


myBATH™ SERIES

MyBath 12L
(B2000-12)



MyBath 4L
(B2000-4)

Tube racks included with myBath 4L!

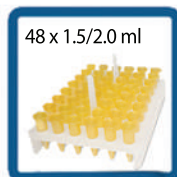

48 x 1.5/2.0 ml



41 x 15ml



15 x 50ml

- Removable, hinged lid*
- Quik-CAL™ user calibration function
- Seamless, stainless steel chamber
- Digital accuracy and convenience

The myBath series is designed to guarantee superior accuracy and performance in a compact footprint. With four sizes to choose from, there is a myBath to fit every lab.

The 4, 8 and 12L models include easy to use controls allow the user to digitally select and monitor settings. Adjustable in 0.1°C increments, the temperature is shown on a large, 3-digit display. The proprietary Quik-CAL™ feature allows for simple, in-lab calibration.* High efficiency insulation maintains temperature stability and uniformity within the bath, while keeping the housing cool to the touch.

All myBath models feature seamless, stainless steel chambers for easy cleaning and long life. The chambers resist rust and have no corners to collect dirt or contamination. The larger 8L and 12L units have drains for hassle-free water changes.

The 4L, 8L, and 12L models include a see-through cover that can be removed, or flipped back on a rear hinge to conserve bench space. The lid's gabled design will divert condensation to the front and back of the bath, away from samples.

Racks are available to accommodate 1.5/2.0ml, 15ml and 50 ml tubes within the 4L, 8L and 12L baths. The myBath 4L is supplied with all three racks.

Technical Data and Ordering Information:

*2L models use traditional calibration and include a flat, removable lid.

	2 Liter	4 Liter	8 Liter	12 Liter
Catalog Number	B2000-2*	B2000-4*	B2000-8*	B2000-12*
Temperature Range	Ambient + 5 to 100°C			
Temp. Accuracy	± 0.2°C at 37°C		±0.5°C at 37°C	
Temp. Increments	1.0°C	0.1°C		
Temp. Uniformity	± 0.2°C at 37°C			
Exterior (WxDxH)	6.5x8x6 in.	7.7 x 13.9x9 in.	11.2 x 17.5x 12.4 in.	11.2x 17.5 x 16.2 in.
Exterior (cm)	16.5x 20.3 x 15.2 cm.	19.5 x 35.3 x 22.9 cm.	28.4x 44.5 x 31.5 cm.	28.4 x 44.5 x 41.1 cm.
Chamber (WxDxH)	5.5 x 6 x 4 in.	5.3 x 9.4 x 5.5 in.	9.5 x 12.8 x 5.9 in.	9.5 x 12.8 x 7.9 in.
Chamber (cm)	13.9 x 15.2 x 10.2 cm.	13.5 x 23.9 x 13.9 cm.	24.1 x 32.5 x 14.9 cm.	24.1 x 32.5 x 20.1 cm.
Weight	4.5 lbs / 2.4 kg	15 lbs / 6.8 kg	17 lbs / 7.7 kg	17 lbs / 7.7 kg
Electrical	115 or 230V, 50-60Hz, 250W-700W			
* 115V with US plug. For 230V (EU plug), please add (-E)				

* 115V with US plug. For 230V (EU plug), please add (-E)

Accessory Tube Racks:

Cat. Number	Tube Rack Desc.	4L Rack Capacity	8/12L Rack Capacity
B2000-4-T5	40 x 0.5ml tubes	1	2
B2000-4-T1520	48 x 1.5/2.0ml tubes	1 (Included)	2
B2000-4-T150	41 x 15ml tubes	1 (Included)	0
B2000-4-T500	15 x 50ml tubes	1 (Included)	0
B2000-8-T150	76 x 15ml tubes	0	1
B2000-8-T500	30 x 50ml tubes	0	1



MyBath 2L
(B2000-2)